

# ENDO MOTOR

CE 0197

Micro-motor for root canal preparation with integrated apex locator

## iRoot pro



➤ **APEX LOCATOR** + **ENDO MOTOR**


### TWO WORKING MODES



type	portable
memory	10memory setting + quick start memory
rotation	clockwise/anti-clockwise/reciprocating
torque control	auto reverse/auto stop
battery	rechargeable/replaceable
charging	charging base/Micro USB
display	OLED/Blue Teeth
available contra angle	4:1 10:1 16:1 20:1
apex locator control motor working mode	auto start/auto stop/auto reverse

The most convenient root canal treatment system , complete root canal length measurement and root canal preparation on the handpiece, you just need to connect a lip hook with the motor to make the canal length measurement and preparation at the same time

Apex locator control motor working mode: Auto start/Auto stop/Auto reverse

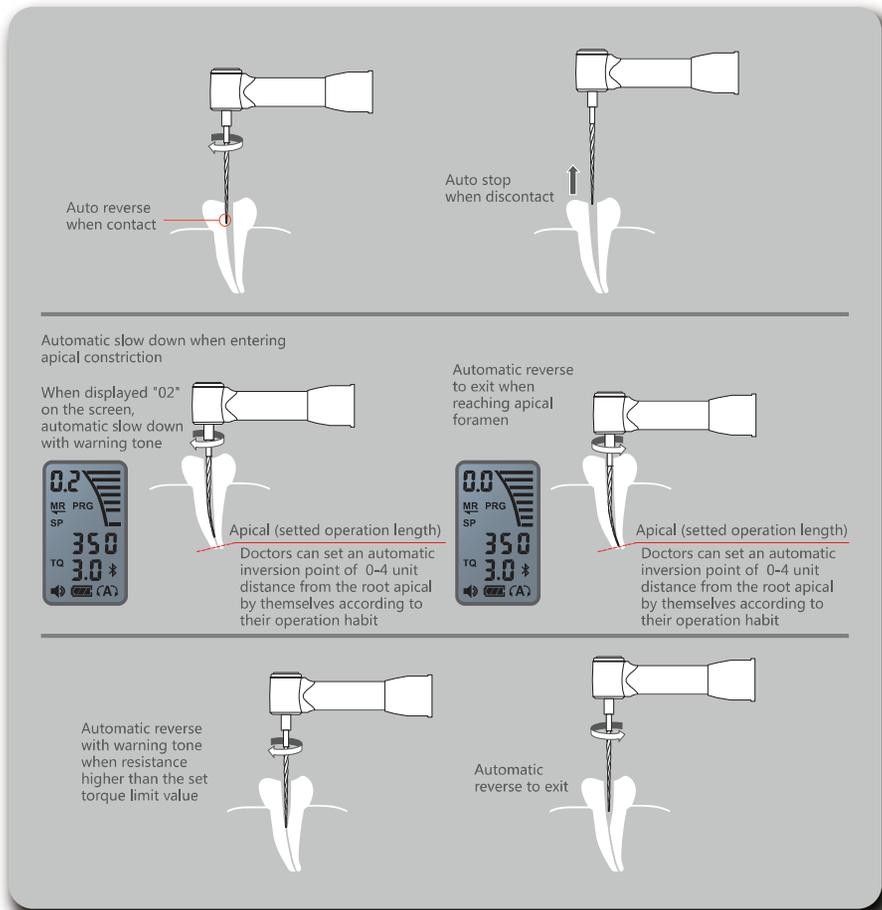

Torque and Speed

The product can be set up with ten memory modes. Auto reverse mode and reverse alarm function make the treatment process much safer. Users can set the corresponding speed and torque according to the parameters of different files at home and abroad. Users can get different speed by matching the contra angle with different reduction gear ratio. The product is applicable to 4:1, 10:1, 16:1 or 20:1 contra angle. For the Reciprocating Rotation, there are clockwise cutting and anti-clockwise cutting. The speed of the reciprocating rotation is adjustable.

The most convenient root canal treatment system, complete root canal length measurement and root canal preparation on the handpiece, you just need to connect a lip hook with the motor to make the canal length measurement and preparation at the same time



We are committed to producing safe and reliable medical equipment to let doctors use more handy, so we recommend doctors to use the apex locator control motor working mode